

REPORT OF ANALYSIS


For: (35432) OPEN FARM -
Dog-OF Ancient Grains Pork
21622FPGD Open Farms Certified Humane Pork Grain

Analysis	Level Found		Reporting			Analyst- Date	Verified- Date
	As Received	Units	Limit	Method			
Sample ID: 21622FPGD-MICRO	Lab Number: 13874993	Date Sampled: 2022-08-05					
Salmonella	negative	org/375g	1	RapidChek/AOAC RI 030301; AFNOR SDI 34/01-04/10	Hct6-2022/08/08	jzh4-2022/08/08	
Sample ID: 21622FPGD	Lab Number: 13874994	Date Sampled: 2022-08-05					
E. coli (generic)	n.d.	cfu/g	10	AOAC 2018.13	Jcp2-2022/08/07	jzh4-2022/08/07	
Sample ID: 21622FPGD	Lab Number: 13874995	Date Sampled: 2022-08-05					
Listeria	negative	org/125g	1	RapidChek/AOAC RI 020401	Hct6-2022/08/08	jzh4-2022/08/08	

This report was reissued on 2022-08-25 10:09:01 by jbw4 for the following reason:
combine reports.

All results are reported on an AS RECEIVED basis., n.d. = not detected , cfu = colony forming unit

For questions please contact:

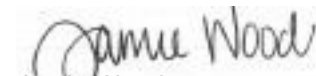

Jamie Wood
Account Manager
jwood@midwestlabs.com (402)590-2964

The result(s) issued on this report only reflect the analysis of the sample(s) submitted.

Our reports and letters are for the exclusive and confidential use of our clients and may not be reproduced in whole or in part, nor may any reference be made to the work, the results, or the company in any advertising, news release, or other public announcements without obtaining our prior written authorization.

REPORT OF ANALYSIS

For: (35432) OPEN FARM -

Dog-OF Ancient Grains Pork

21622FPGD Open Farms Certified Humane Pork Grain

Detailed Method Description(s)

Salmonella - Lateral Flow

Samples are analyzed following MWL MI 195 which is based on the RapidChek Select Salmonella User Guide. A representative sample is obtained and combined with a selective media and allowed to incubate. After incubation, a test strip is used for Salmonella determination. Results are reported as negative or presumptive positive.

E. coli and Total Coliform using 3M Petrifilm

Sample analysis follows MWL MI 292 which is based on AOAC 2018.13. A representative sample is obtained and added to phosphate buffer. Aliquots of the sample are withdrawn and placed on Petrifilm plates. The plates are incubated for 18 to 24 hours. After incubation, the plates are counted to determine the number of generic E. coli and total coliforms present. The color of the colony and the presence of gas differentiate a generic coliform from E. coli. The levels are reported as colony forming units (cfu).

Listeria Lateral Flow

Samples are analyzed following MWL MI 194 which is based on the RapidChek Listeria User Guide. A representative sample is obtained and combined with a selective growth media. It is incubated for 40-48 hours. After incubation, an aliquot is heated for 10 minutes, and a test strip for Listeria detection is used. Results are reported as negative or presumptive positive. This procedure does not speciate Listeria.

The result(s) issued on this report only reflect the analysis of the sample(s) submitted.

Our reports and letters are for the exclusive and confidential use of our clients and may not be reproduced in whole or in part, nor may any reference be made to the work, the results, or the company in any advertising, news release, or other public announcements without obtaining our prior written authorization.